

Product Parameter

Model	Light source	Illuminance	Size of illuminance field	Color temperature	Color rendering index (CRI)	
Whale	4 light-emitting diodes (LED)	Adjustable from < 8,000 to > 25,000 lux according to ISO 9680	80mm×160mm@700m	5,000K-5500K	Ra≥90, Rf≥90	
	On/Off control		Dimming control		Power input	
	Non-touch sensor		Non-touch sensor		AC 12 or AC 24V (Preset before delivery)	
	Power consumption		Degree of protection against ingress of water		Mode of operation	
	< 15W		Ordinary equipment (IPX 0 - without protection against ingress of water)		Continuous operation	
	Transport and storage conditions			Operating conditions		
	Temperature		-40℃~+70℃		Ambient temperature	
	Relative humidity		10%-95%		Relative humidity	
					30%-85%	
	Illuminance diagram			Dimension of light head, movement and mounting		

ADS
ADS Dental System

Whale



Post mount Whale LED Dental Light 30"/40"

Part number:

A0602600 30" Post mount
A0602602 40" Ortho Post mount



Ceiling mount Whale LED Dental Light

Part number:

A0601601 8' Ceiling mount
A0601602 9' Ceiling mount
A0601603 10' Ceiling mount

2020 version1

Technical support: 800 488 9708

Add: 1590 S Milliken Ave., unit A, Ontario, CA 91761

Tel: (626) 620 0456

www.adsdental.com www.adsequip.com

sales@adsequip.com





Whale

LED Dental Operating Light




Perfect experience of the unique features



Comfortable vision

Comfortable, tireless, clear view under NiceDay's Color and Intensity uniform light pattern, unique feather like soft boundary take away the tiresome sharp change of light.

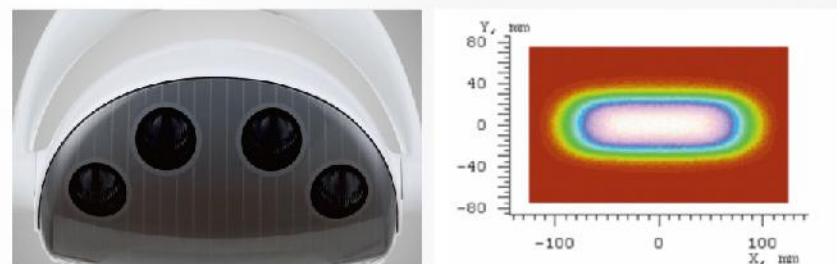

Part number: A121709



Easy to clean

Hygiene, Hygiene and hygiene! Niceday take all the opportunities available for better performance of hygiene. Non-touch sensor for ON/OFF and dimming, Press button release handle for cleaning. Enclosed simple over all shape, as less seams as possible, but still looks so great that will improve the image of your practice!



Smooth & accurate movement

Confident adjustment of the light beam, the all 3-Axes moving hinge is designed to duplicate the accurate and comfortable feeling of the fine turning of high end Camera's lens. Each movement become an Enjoyment!



Long service life

High quality LED from leading international company, driven by state of art IC, Aluminum casting shell as main heat sink, all these guarantee the 30,000 Hours expected service life of the LED light source.

